

PRACTICAL INNOVATION



AUTOMATED PLATE STACKER

Product ID – ZSPS-01



www.ziath.com

System Specification

- Automated SBS format plate stacking system
- All SBS format micro well plates
- Cycle time of 7 seconds present 7 seconds return
- Ability to handle tip racks; contact Ziath to confirm ability of rack types
- Plate heights from 9 to 48mm high
- Handles lidded and non lidded plates (requires optional easy hotel change)
- Ability to de-lid plates and retain lid for replacement (requires optional easy hotel change)
- Holds up to 72 plates
- Supplied with two hotel units: one pre-, one post-cycled plates
- Narrow footprint to minimise working area
- Twin systems can operate to increase in and out feeds to operational units
- Instant interface with Ziath/KBiosystems equipment
- Fully integrated with the following equipment
 - KBiosystems Chameleon / Ziath Automated Sealer
 - Thermo ALPS 300 / ALPS-3000
 - Porvair Tri-Seal
 - Beckman DTX and ParaDIGM
 - Nexus Xpeel
 - BMG Plate readers
 - Ziath A6 2D Scanner
- Easy operation user control software
- Adjustable plate grippers to suit plate changes
- Available in both left and right hand feed
- Offers a restack option if plate arrangement is important
- Integration friendly to other plate based systems (full ASCII command set available)


System Footprint

Maximum Arm Travel

Standard Left and Right Arms: 280mm

Maximum Arm Travel

Extended Left and Right Arms: 342mm

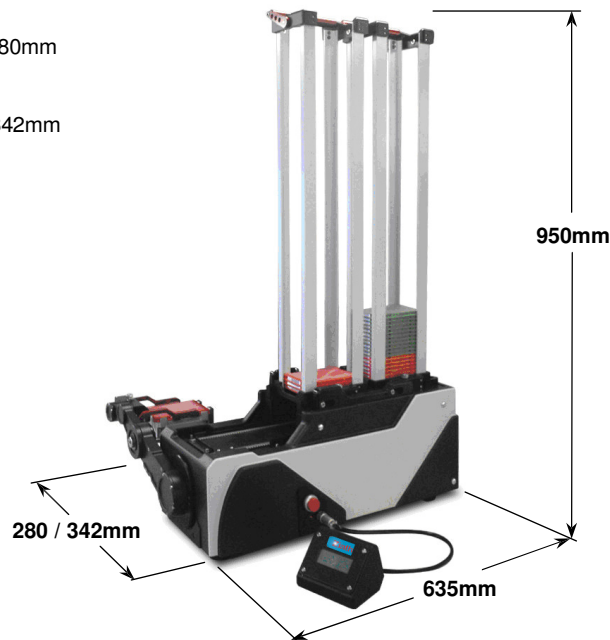

For more information please contact Ziath Ltd.

www.ziath.com • sales@ziath.com

T +44 (0) 1223 968 125 • F +44 (0) 1223 771 615